


# VariAx<sup>®</sup> 2 Plating System



# VariAx 2




## Plating System

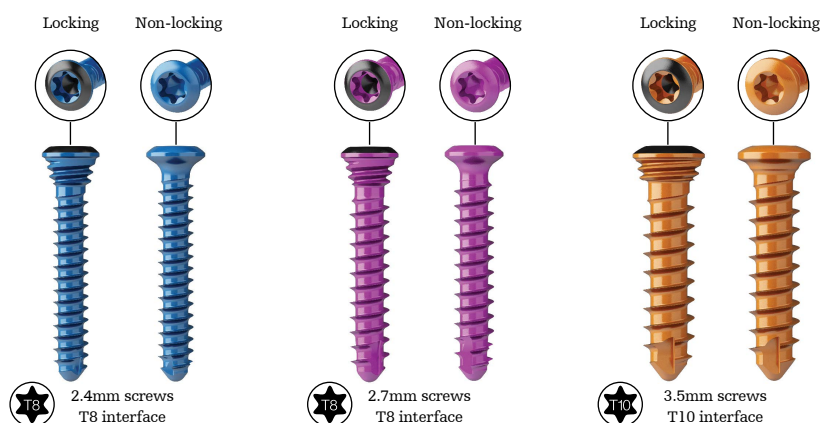

### Comprehensive locking and non-locking system

Our new color-coded polyaxial screws provide a range of length and a new color code for easy identification for both diameter of the screw and between locking and non-locking bone screws. Each VariAx 2 foot plate screw hold incorporates patented SMARTLock (licensed from Prof. Wolter) technology. Each hold also accepts non-locking screws. With a single plate comes more options. The genius of SMARTLock is its simplicity.

The VariAx 2 screws and instruments follow a comprehensive color-coding scheme whereby the screw anodization color matches the corresponding instrument's color. This allows for easy handling in the OR environment. Overdrills and drill guides for interfrag-mentary lagging technique are available for all screw sizes. The drill bits, lagging drill guides and screws are also color coded to facilitate ease of use. Drill guides and screw-drivers are also identified with a T10 or T8 symbol for quick reference.

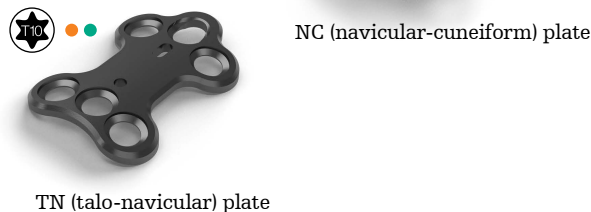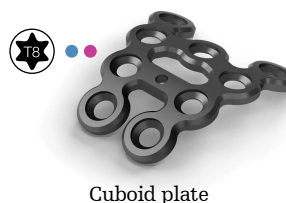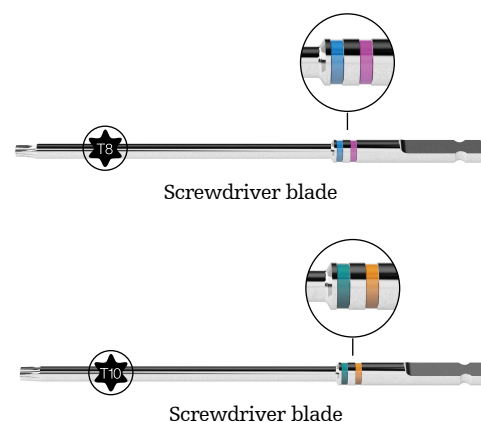

TN (talo-navicular) plate

### Foot & Ankle

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate a Stryker product. A surgeon must always refer to the package insert, product label and/or instructions for use, including the instructions for Cleaning and Sterilization (if applicable), before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: SMARTLock, Stryker, VariAx. All other trademarks are trademarks of their respective owners or holders.

VAX-BR-7 Rev-1, 08-2016  
Copyright © 2016 Stryker



**Manufacturer:**  
Stryker GmbH  
Bohnackerweg 1  
2545 Selzach  
Switzerland

stryker.com